


## Emulator - Type 6

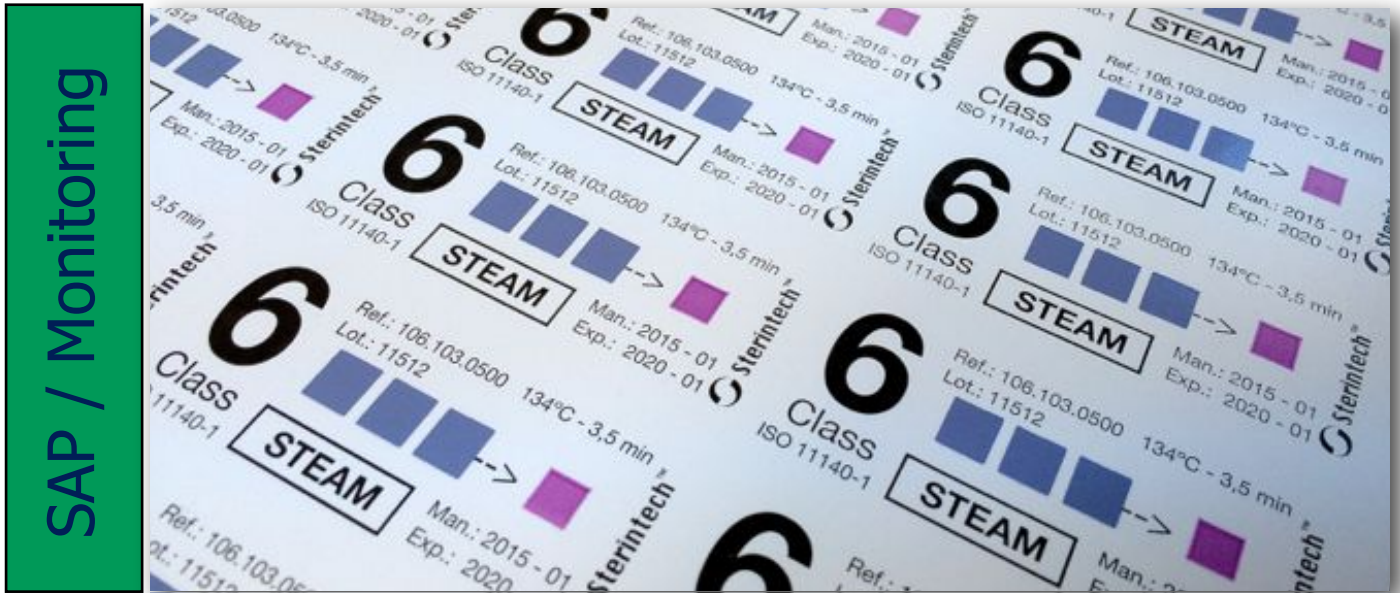

### The Emulator with a big plus...

Emulators can be used in addition to the batch control Helix which will allow you to release load. Especially when OR staff wants to be re-affirmed about the sterilization conditions.

REF: 106.103.0500SA	Emulator 134°C-3,5 min/121°C-15 min Includes Self-Adhesive backing	Unit: 500 pieces
REF: 106.103.0500	Emulator 134°C-3,5 min/121°C-15 min	Unit: 500 pieces

The Sterintech™ Emulators are offering you more information about the steam sterilization process by means of the very accurate indicator fields of the Emulator. These Emulators can detect many disturbances of the sterilization parameters inside packs or containers.

Characteristics:

- ◆ Shelf life : 5 years after production
- ◆ Steam emulator conform to ISO 11140 - Type 6
- ◆ Clear color transition from Blue into Pink
- ◆ Chemical emulator covered with lacquer
- ◆ Non-Toxic / Lead Free




Majac Medical Products Pty Ltd

Ph: +617 3265 6355  
Fax: +617 3865 2729  
sales@majacmedical.com.au  
www.majacmedical.com.au



FS 540017

Also Available From:

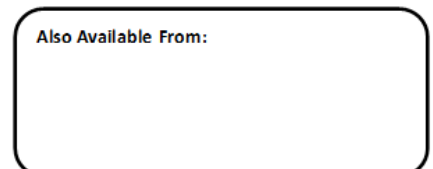

## Sequential Emulator - Class 6



### The Emulator with a big plus...

Integrators can be used in addition to the batch control Helix which will allow you to release load. Especially when OR staff wants to be re-affirmed about the sterilization conditions.

REF: 106.003.0250	Sequential Emulator 134°C-3,5 min	Unit: 250 pieces
REF: 106.007.0250	Sequential Emulator 134°C-7,0 min	Unit: 250 pieces

The Sterintech™ Emulators are offering you more information about the steam sterilization process than those with only one emulator printed. Due to that the three emulators are calibrated to different time/temperature combination it is providing information about for instance over-heating.

Characteristics:

- ◆ Shelf life : 5 years after production
- ◆ Packed in an easy dispensing carton
- ◆ Steam emulator conform to ISO 11140 - class 6
- ◆ Clear color transition into black of the emulator
- ◆ Chemical emulator covered with lacquer
- ◆ Non-Toxic / Lead Free

