
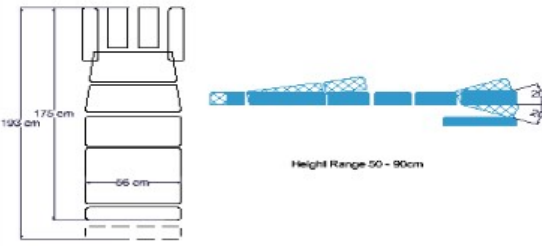


Elevation Chiropractic Table with Drops

	<p>The ABCO Elevation Chiropractic Adjusting Table with Drops is designed and manufactured in Australia.</p> <p>With Elevation [height adjustment] function handling up to 250kg, this is a versatile table providing flexibility for your practice</p> <p>The twin column structure:</p> <ul style="list-style-type: none"> • Provides a much more rigid working surface without 'end-bounce' • Eliminates lengthwise movement during height adjustment
	<ul style="list-style-type: none"> • Improves the safety of the couch as there are far fewer points where stray fingers or limbs may be caught • Enhances treatment by moving quickly and quietly as required <p>You will love our crisp, heavy duty, fully adjustable drop pieces. Four drops are included with hand cocking as standard. Foot cocking is available as a no-charge option</p> <p>Product Code: CIVD.05642.CDDDDF.50NP0</p>

Detailed Standard Specification

Table Sections	Cervical Section	320mm long with drop. Adjusts from 15° to -25°
	Dorsal Section	240mm long with drop
	Lumbar Section	210mm long with drop
	Pelvic Section	240mm with drop 560mm long. Adjusts to 10
	Leg Extension	100mm long; extends 100mm
	Overall Length	1900mm with leg section extended
	Upholstery angle adjustment	Gas strut assisted
Drops	Drop Height	16mm
	Drop Force	0 - 100% adjustable from either side of table
	Mechanism	Linear bearing / ball bearing
	Cocking Mechanism	Roller bearing from either side of table. Hand operated
Upholstery	Fabric	Leather - New Zealand Pigmented heavy commercial grade cattle hide.
	Foam	80mm very high-density foam with 12mm soft foam over
	Standard Upholstery Colours	Black. Other colours on request
Frame	Materials	Heavy steel plate and tube
	Colour	Light grey
	Wheel System	No wheels
Height adjustment	Lift Mechanism	Twin column for maximum rigidity in use.
	Lift Operation	Linak brand electric columns
Electrical requirements		Medically approved according to EN60601-1/ AS3200
	Supply	230 VAC +/-10%, 50Hz max 1.5amp.
	Operation	24 VDC max 70VA
	Control	Foot control
Maximum lift		> 300 kg
Table weight		130 kg
Shipping Dimensions	Standard product	1800 x 600 x 650mm / 145kg
		Shrink-wrapped on a wooden pallet with bubble wrap
ARTG No:		123573
Warranty		5 years, parts only