Quartz (AMC 2530)

A general scanning couch with a cut-away section, allowing the transducer to be positioned more easily for some cardiac views.

Electric adjustments of height and back rest.

Designed to reduce OH&S risks and minimise manual handling.



warranty

on electrical components



warranty on frame and workmanship







*Pictured central wheel locking not included as standard. Available as an optional customisation.

Features

- > Adjustments of height (hi/lo) and back rest made with high quality electric actuators, reducing OH&S
- > Supplied with foot control switch for hands free operation.
- > Permanent opening for cardiac scanning.
- > Carefully considered base design allows easy side access for practitioners and patients.
- > Access under couch for patient lifting devices.
- > Low minimum height the gives patients the ability to easily get on and off the couch.
- > 125mm lockable wheels allow easy movement of the couch.

Basic specifications

Width: 700mm Length: 1970mm Back rest tilt: 75°

Height adjustment: 550mm -1000mm

Lifting capacity: 200kg

Wheels: 125mm lockable castors

Optional customisations

- > Hand switch control
- > Battery pack
- > Head cushion
- > Paper roll holder
- > Wide variety of fabric colours available
- > Padded swivel arm rest set
- > Standard arm rest set
- > Fold down side rails
- > Central wheel locking



Majac Medical Products Pty Ltd


Ph: +617 3265 6355 Fax: +617 3865 2729 sales@majacmedical.com.au www.majacmedical.com.au




Also Available From:

Quartz (AMC 2530)

Side view

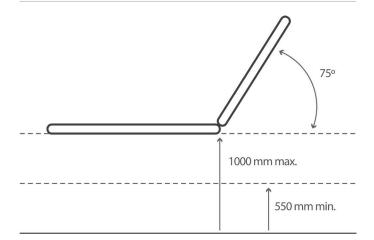

Top view

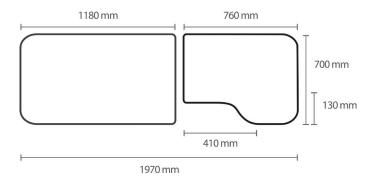

Optional customisations



Head cushion



Central wheel locking



Fold down side rails



Paper roll holder

Detailed specifications

Head section length: 760mm

Middle section length: N/A Foot section length: 1180mm

Vinyl specifications:

UV resistant Flame resistant Easy to clean

Reduces growth of bacteria Complies with AS 1530.3 standards

Variety of colours available

Foam specifications: 50mm foam

Frame material and colour:

Powder coated white frame

Australian steel

Wheel system:

Individual 125mm brake castors made by Tente

Drive actuators:

Back rest: Linak electric actuator system Hi/lo: Linak electric actuator system

Electrical specifications: 230V 50-60 Hz

111/10. Illiak electric actuator system

Electrical standards: IPX4 rating

Lift capacity: 200kg

Shipping weight: 70 - 85 kg

Shipping dimensions:

2130mm x 820mm x 700mm

Freight charges:

Not included in sale price.

Contact our sales representatives for a quote.

ARTG number: 208309

Warranty:

Lifetime on frame and workmanship

5 years on electrical components



Majac Medical Products Pty Ltd

Ph: +617 3265 6355 Fax: +617 3865 2729 sales@majacmedical.com.au www.majacmedical.com.au



Also Available From: