


Batch Control Helix Test (VH₂O₂)

SAP / Monitoring

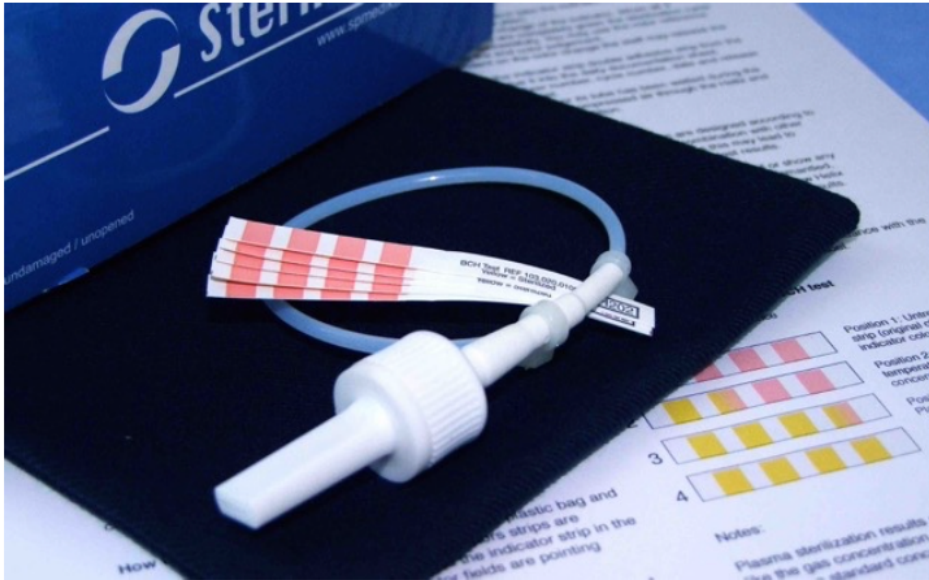

Checks your cycles... on all parameters

REF: 103.020.0100

For plasma cycles

Unit: 100 indicators

Helix devices and their indicators have proven that they detect failures during each cycle. From leakages, too less plasma, overloading, time and vacuum issues and forms a simple aid to release your load securely.

The PCD is consisting of a 30 cm long tube with an internal diameter of 2mm which is the maximum internal lumen which can be sterilized effectively by plasma sterilizers (single open end)

Characteristics:

- ◆ Conform to ISO 11140 - 1: 2014 Type 2
- ◆ Clear transition of color into yellow
- ◆ Easy to interpret
- ◆ Easy to document
- ◆ Non-Toxic / Lead Free
- ◆ Can be fully recycled
- ◆ PCD guaranteed to min. 100 cycles



Sterintech™



MAJAC
HEALTHCARE INNOVATION

Marketed by:
Majac Medical Products Pty. Ltd.
Tel: 1300 138 578 • Int: +61 7 3265 6355
Email: sales@majacmedical.com.au
www.majacmedical.com.au



FS 540017

"Let's Get it Right"