

# Bovie®

## SPECIALIST|PRO

*Electrosurgical Generator*



The Bovie® Specialist Pro is the electrosurgery solution for in demand physician offices and specialty facilities. Designed to operate like an operating room generator in a compact lower watt package. The multipurpose unit is ideal for monopolar and bipolar office surgical procedures.

The Specialist|PRO offers 5 common energy outputs controlled by the new sleek, push-button design, with up to 120 watts of power. The modalities include Cut, Blend, Coagulation, Fulguration and Standard Bipolar.

***Efficiency, Versatility, and Safety ...***  
*Designed benefits for you and your patients*

Industry-standard connections allow for use of a variety of monopolar and bipolar accessories providing energy versatility across a variety of specialties such as OBGYN, Urology, Plastic Surgery, Dermatology, ER and more.

Bovie NEM™ neutral electrode monitoring increases patient safety by ensuring proper pad connectivity during procedures.

The Bovie Specialist|PRO is backed by the unmatched **4-year manufacturer warranty**, offered only by Bovie.



Visit [BovieMedical.com/specialistpro](http://BovieMedical.com/specialistpro) for more information on the Bovie Specialist|PRO.

# BOVIE® SPECIALIST|PRO ELECTROSURGICAL GENERATOR



## STANDARD EQUIPMENT

- A1250S electrosurgical generator,
- ESP1 sterile disposable pencil
- ES01 blade, ES02 needle, ES20 ball (1 each)
- sterile disposable electrodes
- ESRE (5) split disposable return electrode
- A1252C (2) reusable grounding cords
- 110 VAC Hospital-grade power cord\*
- 220 VAC Hospital-grade power cord\*
- User's guide (MC-55-237-001)\*

\* Not Pictured



ESP1



ES01, ES02, ES20



ESRE



A1252C

## COMPLEMENTARY PRODUCTS †



SE02  
Smoke Evacuator †



BV-1253B  
Monopolar Footswitch †



A1255A  
Footswitch Adapter †




Disposable LLETZ/LEEP  
Loop and Ball Electrodes †



ESREC - Split Disposable  
Return Electrode with cord †



Disposable SuperCut™  
Needles Electrodes †

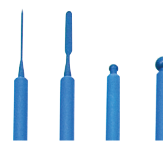

Disposable Resistick II™  
Electrodes †



Disposable Dermal Loops †

† Sold Separately

## PRODUCT SPECIFICATIONS

### DIMENSIONS

Height: 15.2cm (6 in.)  
Width: 26cm (10.25 in.)  
Depth: 30.5cm (12 in.)  
Weight: < 4Kg (< 9 lbs.)

### OUTPUT POWER

Cut: 120 watts  
Blend: 90 watts (60% Cut, 40% Coag)  
Coagulation: 80 watts  
Fulguration: 40 watts  
Bipolar: 30 watts

### OUTPUT FREQUENCY

350 – 600 kHz

### LINE VOLTAGE

A1250S: 100-240 VAC

### LINE FREQUENCY

50 – 60 Hz

### WARRANTY

Four Years



Marketed by:  
**Majac Medical Products Pty Ltd**  
Tel: 1300 138 578 \* Int: +617 3265 6355  
Email: sales@majacmedical.com.au  
www.majacmedical.com.au



FS 540017

Also available from: